

# Im.PAc.T


## meticule SBRT SYSTEM

Model No. : ROOSBRT30-00




Panacea's Im.PAc.T Meticule SBRT System is effortless to use and is indexable on any couch top available at your clinical facility. The system is entirely made up of carbon fiber which increases its durability over the years of its usage. The standard SBRT system consists of rails that are indexable on all the couch tops, extended wing board and adjustable respiratory bridges which aids in the most comfortable positioning for radiation treatment.


To provide uniform compression, Im.PAc.T Meticule SBRT system comes with a respiratory belt which can be customized to the patient's body size.

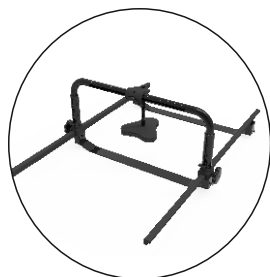


The movable rails can be indexed directly on all the couch tops. The respiratory bridge is adjustable to patient's body size with vertical and transverse travel range up to 28cm and 36cm respectively ensuring comfortable patient positioning.

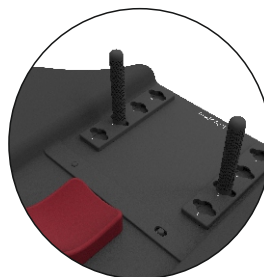

The unique extended wing board can be indexed on the couch top directly and consists of adjustable and secured hand grips to ensure maximum support to the patient.

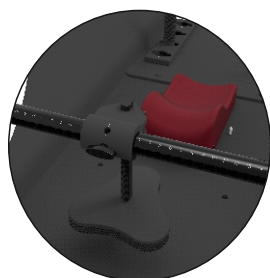

The respiratory bridge is equipped with a soft cushioned compression paddle which is adjustable up to 7cm while compressing the abdominal region without causing pinching and pain to the patient.

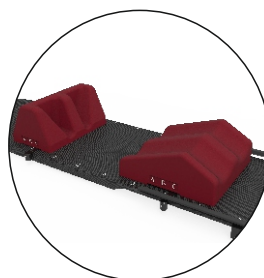

The knee fix and foot fix has three indexing positions to raise or lower the height for maximum leg comfort. It is made of skin safe foam cushions and are easily mounted on couch top/ SBRT Board.